


## Changes in the axis sleeves of caster fork on Agile

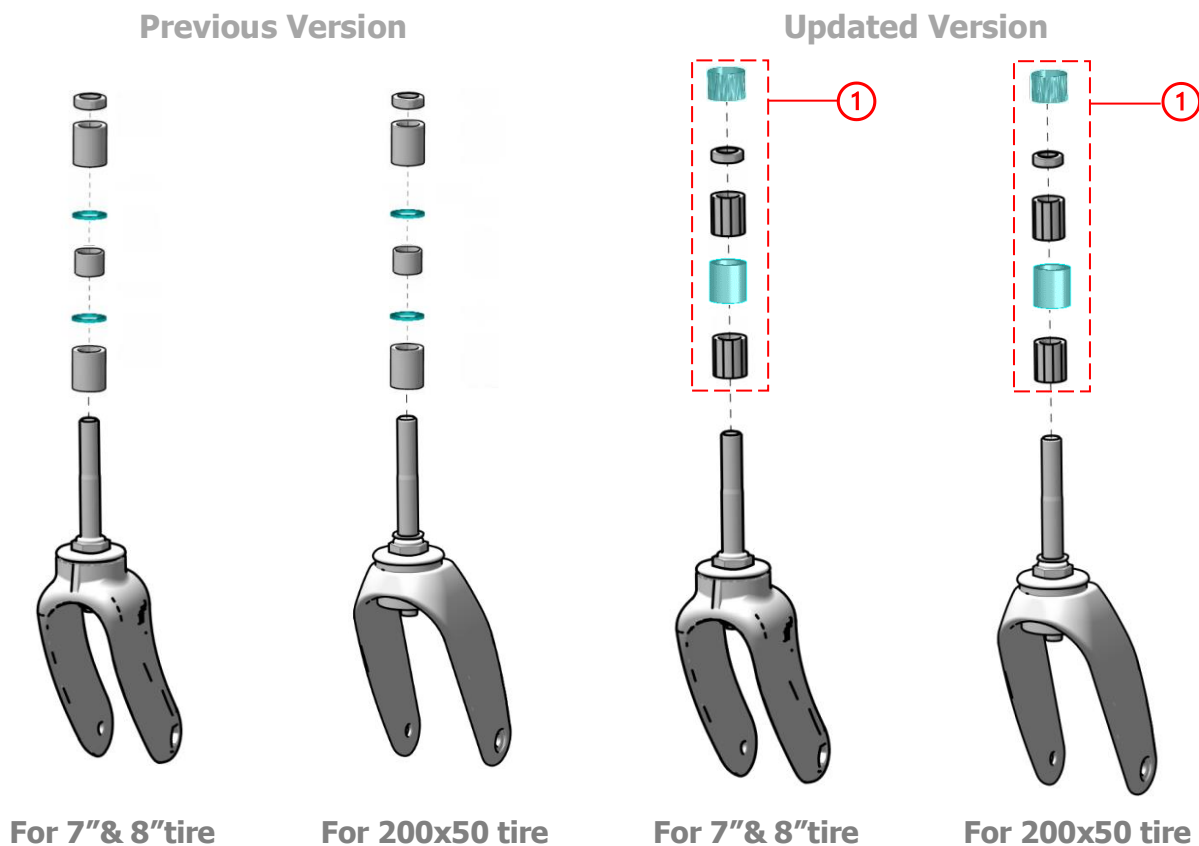

- ✓ To improve the tightness between the caster fork axis and side frame on Agile, Karma has changed the design of the axis sleeves as follows:
  1. Two flat washers have been removed from the axis sleeves kit.
  2. The material of the middle axis sleeve has been changed from aluminum to plastic.
  3. The new anti-slip plastic sleeve has been added on the nut that secures the top of the caster fork axis.
- ✓ From now on, please order the new sleeves kit ① when repairing the front caster fork.
- ✓ The updated version is backward compatible with the existing one.
- ✓ The updated version was implemented by the batch of serial no. HM22-025623, not including HM22-27390~HM22-27404, HM22-27375~HM22-27389 and HM22-27675~HM22-27704.
- ✓ The updated version was implemented on Agile.

### Order Information

Article number	420400000136	FrontCaster Fork Assy 63W*148H*125L
	420400000137	FrontCaster Fork Assy 73W*227H*125L
	① 420400000205	Sleeve kit