

KARMA

Karma Flexx Junior

Specifically designed for functionality
with multiple adjustments to suit
children's needs.



With growth
incorporated



**ALLIED
MEDICAL**
LIMITED

+ FLEXIBLE, ROBUST AND STYLISH

The Karma Flexx Junior is designed with plenty of adjustments to create flexible sitting to fit individual needs. With ellipse-shaped frame design, the Flexx is not only stylish but also robust and durable.

The configurable options of the Flexx Junior include variable seat depths, seat to floor heights, backrest heights and backrest angles, castor fork angles and centre of gravity. It can easily be set up to provide optimal comfort and performance.

The bright colour of key adjustment points are self-explanatory; users can easily operate the following functions independently: swing-away footrests, flip-back and adjustable armrests.

+ FEATURES

- Configurable dimensions
- Height adjustable flip-back armrests
- Swing-away 90 degree footrests with flip-up footplates
- Back angle adjustment (-4° to +16°)
- Padded height and angle adjustable backrest
- Reinforced double crossbrace
- Quick release pneumatic rear wheels
- Push to Lock Brakes with brake extension
- 2 point lap belt
- Heel loops



+ ACCESSORIES

- Swing-away height adjustable anti-tippers
- Tray

+ SPECIFICATIONS



Seat widths
10" (25.5 cm) to
14" (35.5 cm)




Seat depth
11" (28 cm) to
15" (38 cm)




Wheels
5" (front)
20" (rear)



**Seat to floor
height**
14" - 22"
(35.5 - 55.9cm)

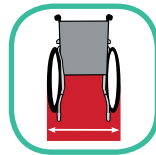

Overall width
18" - 22"
(45.4 - 55.6cm)

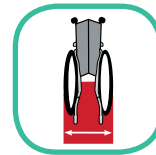

Folded width
12" (30cm)



Back height
14.5" (37cm)



Overall length
35" (89.6cm)



Overall height
32" (81cm)



Total weight
15.7 kg



User weight limit
80 kg

