

Little Wave Flip Accessory Flip Chart

 Ki Mobility



Flip the script on ordinary

Put a flip to your fold

The Little Wave Flip integrates design elements from strollers into the durability of a wheelchair. Just because it is a wheelchair, does not mean it needs to look like one. From its flat folding design to its use of TAPER LOK technology, the Flip is simply the most advanced paediatric tilt in space available today.

Eliminate bland. Illuminate personality

Children are small in size but their personalities can be bigger than life. Let them express themselves with over 500 colour combinations to choose from-no charge!

Designer Spoke Guards

Eight NEW spoke guard designs offer a few more ways to personalize your Flip! View all designs in the "Options" tab above.

Don't struggle with growth. Embrace it

The Flip offers a new paradigm in adjustability. Seat height, seat depth, and seat width have never been so easy to change, with fewer parts than any other chair.

Attendant Wheel Lock

The attendant wheel lock is a hands-free, easy to use stroller-like wheel lock option that makes parents life just a little simpler.

Ergo Folding Back

The Ergo Folding Back is tough, strong and easy to fold. With the flip of a lever you are able to activate the Ergo Folding Back for an effortless fold.

High Pivot Folding Back

By raising pivot point on the folding back, we've created the possibilities to leave more seating on the chair while folded.





Frame

Transport Tie Downs



Conforms to WC-19 Standard.

Offset Hanger



Available only with Rigid Base

Tilt Mechanism



Available as Hand or Foot.

Backrest

Stroller Handle Back Post



Available in Short (20") or Tall (24"), the back cane is vertical to the frame with a 90degree bend at the top for a push handle. Folding

Fixed Height Back w/ Adj Handle



1" diameter to accommodate a multitude of aftermarket seating.

Adj Height Back W Adj Handle



Available in Short and Tall height ranges, this back post is height adjustable and comes with an adjustable handle. Folding

Stroller Handle Extension



A plug-in extension to assist in the proper placement of the grips for the caretaker.

Seating Dynamics Rocker Back



Shock absorbing back frame reduces stresses on the frame from high tone riders.

Not available with Flip Up Armrests. Available only with Non-Folding Back.

Adjustable Height Rigidizer Bar



Allows you to reposition the rigidizer bar.

Footrest

Extension Tube Hanger



The most intuitive hanger on the market allows you to attach and detach from any angle or position. Available in 60°, 70°, and 80° versions. Short, Medium, Long.

Front Mount Hanger



With 3.5" knee-to-heel capability in .5" increments. Available in 70°, 80° and 90° versions.

Elevating Leg Rests



Simple to use, the elevated leg rests raise up on command and lower easily by using the lever on the side. Also available in extended length version.

Pro Elevating Leg Rests



Pro ELR has an adjustable pivot location to promote proper knee placement and an adjustable calf pad for personalized fit from client to client.

Composite



Made from a fiberglass filled nylon composite, these footplates are lightweight and tough, making them hard to bend or break.

Composite Angle Adjustable



The Composite Angle Adjustable footplate makes adjustments fast and easy. Includes adjustments of the angle as well as fore and aft positioning of the footplate.

Aluminum Angle Adjustable



Infinite adjustability for more accurate placement of the foot. Aluminum plates have slots for ankle straps and multiple attaching holes.

One Piece Flip Up Angle Adj



Solid, angle adjustable, aluminum footplate provides an excellent surface for your feet and increases rigidity at the front.

Heel Loops



Cradles the foot and helps keep it on the footplate. Not available with One Piece Flip Up Angle Adjustable Footplate.

Stump Support



Just as easy to operate as our patented swing away hangers using the same swing in/out design. An incredible amount of adjustment within the system allows you to achieve the exact position and angle your client needs. Available in Smooth or Air Mesh.

Velcro Adjustable Calfstrap



Sitting behind your calf it is adjustable in length and supports your lower leg.

Padded Velcro Calfstrap



This calfstrap supports a soft and durable 4"-wide neoprene pad.

Bodypoint Aeromesh Calfstrap



Provides comfortable support, keeping feet safely on the footplates. Retains its shape for the best pressure distribution and durability.

4 Way Hanger Release - Inline



Mount in 4 different positions with ONE bolt! 4 - Way - In Line position provides best optimization for caregivers.

4 Way Hanger Release - Out



Mount in 4 different positions with ONE bolt! 4 - Way - Out provides best optimization for clients.

Casters

.75" Lighted Rollerblade



Similar to the standard microcaster. Add some pizzazz to your chair!
Not available with 6" and 8" caster sizes.

1" Poly



This outstanding design uses molded polyurethane to create a great tyre. Measuring 1" in width, this caster is available in 4", 5", 6", 7" and 8".

1" Poly Aluminum



Uses a 3-spoke aluminum hub for added strength and styling.
Available only with 4" and 5" Casters.

1" Pneumatic



A low cost way to get a great and comfortable ride. You just need to add air from time to time.
Available only with 6" Casters.

1.5" Poly



This outstanding design uses molded polyurethane to create a great tyre. Measuring 1.5" in width, this caster is available in 4", 5" and 6".

1.5" Soft Roll Aluminum



3-Spoke ultralight aluminum hubs, highly resistant, with improved mobility. Available in 5" and 6".

2" Poly



This tyre is quite firm with a subtle radius so that it is easy to push on hard surfaces and gives you great flotation on soft or uneven ground.
Available only with 6" Caster Size.

2" Pneumatic W Foam Insert



A low cost way to get a great and comfortable ride. These casters have a foam insert so no air is needed.
Available only with 8" Casters.

Quick Release Axle



Standard axle with a quick release button for easy wheel removal.

Standard Axle Plate



Self Propel Axle Plate



Wheels

Mag



Our 5-spoke composite wheel is 20% lighter than the competition and maintenance free. 12", 16" 20" 22' and 24".

Colored Mag



Our 5-spoke composite wheel, now available with color options in 20" and 22" sizes, is 20% lighter than the competition and maintenance free. Colors available are Pink, Blue, White and Fluorescent Red.

16" Mag Wheel



The new 16x2" mag rear wheels make it easier to handle everyday obstacles such as curbs and inclines. Available in pneumatic, pneumatic with airless insert, full poly tire, or low profile poly.

12" Mag Wheel



Spoke



The Ki standard spoke wheel uses a very high quality machined lightweight aluminum hub and 30 stainless steel straight pull spokes for increased stiffness. The radial spoking increases the stiffness of the wheel. Not available with 12" or 16" Rear Wheel Sizes.

Full Profile Poly 1 3/8"



Corded to prevent it from rolling off the wheel. Provides a softer ride than other low cost "over-molded" tyres.

Low Profile Poly 1 3/8"



Pneumatic 1 3/8"



This standard air tyre is more forgiving when going over obstacles and off curbs or stairs.

Pneumatic W/ Airless Insert 1 3/8"



Low density foam is inserted into a pneumatic tyre. Low maintenance - no air, no flat tyres!

Aluminum Anodized Handrim



Aluminum is lighter than steel and stronger than plastic and won't rust.

Plastic Coated



For people who need assistance with grip, but do not require projections. Dipped in Vinyl, they are coated to create a tacky surface for gripping.

Projection



8-3/8" poly-coated vertical projections for individuals with limited grip strength. Available on 22 and 24" Spoke wheels and 24" Mag wheels with full tires.

Natural Fit LT Regular Grip



Eliminates pushing on the tyre and eases pain in the hands and wrists. Offers a better grip for a better push and more control when braking. The LT is a smaller gripping surface.

Natural Fit Standard Reg Grip



Eliminates pushing on the tyre and eases pain in the hands and wrists. Offers a better grip for a better push and more control when braking.

Natural Fit LT Super Grip



Eliminates pushing on the tyre and eases pain in the hands and wrists. Offers a better grip for a better push and more control when braking. The LT is a smaller gripping surface. The Super Grip uses a higher friction surface for extra grip.

Natural Fit Super Grip



Eliminates pushing on the tyre and eases pain in the hands and wrists. Offers a better grip for a better push and more control when braking. The Super Grip uses a higher friction surface for extra grip.

The Natural-Fit Oval



The Natural-Fit LT Oval



Micro Grip



Micro Grip Non-Slide Tape enhances your grip in dry, wet, and snowy conditions! Made of micro-injected silicone and 3M adhesive. Not available on Plastic Coated or Projection Handrims.

Wheel Locks

Push or Pull to Lock



After engagement the lever is at 45 degrees. Lever is easily reachable.

Push or Pull With Extensions



For people with limited strength, this design requires very little effort to engage. Removable 6" fold-able extension increases leverage.

Attendant Footlock



Drum Brake




Looking for more control while maneuvering on slopes? The drum brake is a true braking system that will give you just that. Available in 16" Poly or Pneumatic

Grade Aids (Hill Holders)



Grade aids prevent your wheel from rolling backwards when you push up a hills.

D Lock's

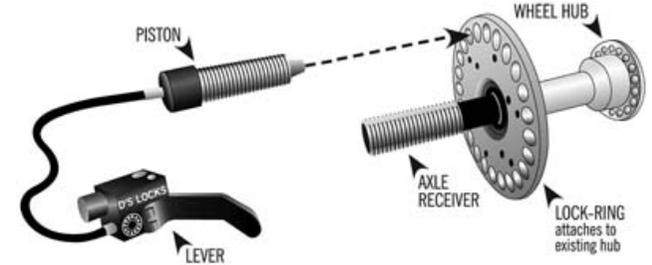

Mounted to your wheelchair and wheel set, D's Wheelchair Locks are a positive wheelchair wheel lock. With D's Wheelchair Locks, all it takes is a simply lever flip to set your wheels in a fully locked position, and they will not turn for anything.

Armrests

Height Adjustable Paediatric Desk



Our single post height adjustable T-arm made for smaller riders. Available in 6.5-9" or Tall 8" - 10.5".

Height Adjustable T-Arm



The single post height adjustable arm is easy to release and adjust in height. Lower profile and lighter than competitive models. Available in 8.5-13.5".

Height Adjustable T-Arm Low



The single post height adjustable arm is easy to release and adjust in height. Lower profile and lighter than competitive models. Available in 7.5-11.5".

Angle Adj Locking Ext Flip Up



This new Angle Adjustable Locking Extendable Flip Up armrest is the ONLY cantilever armrest approved to 400 lbs.! It locks into position by a spring loaded latch while extending up to 4" in 1" increments and 7 angle adjustments in 5° increments from -5 to 25°.

Tubular Flip Up Armrests



Simple Flip up arm engages at the top to stay up for transfers.

Standard Armpads



Padded and textured armrest pad available in both full and desk lengths.

Waterfall Armpads



A generously padded waterfall design providing additional comfort and protection. Available in full and desk lengths. Available only with Angle Adjustable Locking Extendable Flip up Armrest.

Foam Grip



Foam Grip armrest pad available in both full and desk lengths. Available only with Angle Adj Locking Extendable Flip Up Armrests.

Positioning Belts Full range available

1" Auto Buckle



Using a double-pull to ensure the buckle remains in the proper position and attaching hardware for proper placement on the frame our new belts are an industry best value

1" Padded Auto Buckle



1.5" Auto Buckle



Using a double-pull to ensure the buckle remains in the proper position and attaching hardware for proper placement on the frame our new belts are an industry best value.

1.5" Padded Auto Buckle



This belt has a wide 3.5" pad for riders who may need additional padding.

2" Airline Buckle



This belt uses an airline style buckle that can add a degree of security. It uses a wider 2" belt.

BodyPoint Evoflex Pelvic Stabiliser




Stays where you need it. Provides stability and flexibility. Versatile and easy to use. Stiffened end straps keep belt from twisting. Sized for users 15kg to 130kg. Available in Extra Small, Medium and Large.

BodyPoint Belt Mounting Kit



Mounting hardware for the BodyPoint Pelvic Stabiliser.

Seatbelt Clamp

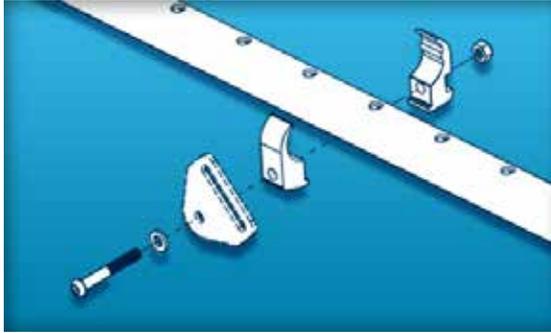

Our new seatbelt clamp allows for 100% infinite adjustment. Ensure proper positioning by achieving the perfect seatbelt angle.

Accessories

Composite Side Guards



Lightweight and attractive, these easily removable side guards help keep road grime away from the rider. Plastic sideguard is height adjustable 14"-20". Pediatric version is height adjustable 10"-14".

Growing Seat Pan



Fixed or removable aluminum seat pan.

Anti Tips



Our design allows the anti-tips to stay connected while making every day obstacles such as curb climbing hassle free.

Tool Kit



Custom designed kit includes 7 simple tools necessary to make any adjustments on your Ki chair, a convenient LED light and carabineer clip for easy attachment.

Removable Underseat Pouch



8.5" x 6.5"



Stylish, functional and accessible backpack for use on chairs down to 9" wide! This smaller size features additional reflective piping details so it's a little more noticeable at night, padded protective pocket, 2 water bottle holders, 5 compartments and adjustable straps.

Backpack



Designed by experts to offer you the most accessible, functional and secure travel storage option. A total of 7 storage compartments. Additional 1.5" depth accommodates a zippered padded storage pocket to hold a laptop or tablet.

Neoprene Impact Guards



Sometimes the chair gets "banged" around. These cool neoprene impact guards take some of that force and help keep your chair a little safer.

Tilt Stops



Spoke Guards



Black or clear acrylic protects fingers from getting caught in the spokes and looks sleek. Available only with Spoke Wheels.

Designer Spoke Guards

