

# KARMA

Karma VIP 515 Tilt-in-space Wheelchair

**Redefining the concept  
of comfort**



**2 Year  
Warranty  
on Frame**

**Folding  
Tilt-in-space  
Wheelchair**



**ALLIED  
MEDICAL  
LIMITED**

## USE TILT-IN-SPACE TO IMPROVE POSITIONING

The VIP-515 is a revolutionary lightweight folding tilt-in-space wheelchair, with a tilt range from 0 to 35 degrees. With the same quality and design features of the well-known Karma 8000 series, it offers a cost effective tilt-in-space option. It is ideal for people in aged care facilities or within the community. With ergonomically designed seating that incorporates some cushioning, it can be used without additional seating products. Its contoured headrest also offers additional support. With a lightweight, folding frame it is easily transportable.

## FEATURES

- Tilt-in-space from 0 to 35 degrees
- Foldable and compact
- Height adjustable and flip-back armrest
- Swingaway legrests that can rotate both inwards and outwards to allow for easy transfers
- Adjustable length legrests
- Ergonomic headrest
- Antitippers
- Heel loops on each legrest



## SPECIFICATIONS



**Seat Width**  
16", 18"



**Seat Depth**  
17"



**Overall Width**  
25", 27"



**Overall Length**  
43"



**Overall Height**  
53"



**Seat to Floor Height**  
20"



**Back Height**  
18"



**Rear Wheel**  
14"



**Front Castor**  
7"




**Armrest Height**  
9"-11"



**Total Weight**  
18.9 kg



**User Weight Limit**  
115 kg

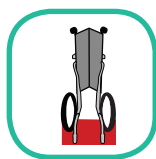

**Folded Width**  
14"



**Tilt Range**  
0 - 35°



**0800 31 61 81**  
[www.alliedmedical.co.nz](http://www.alliedmedical.co.nz)